

# QA3™

## Patient Stretcher

Compact, light weight and maneuverable –  
the complete solution for patient transport



Innovative **Medical** Technology - **Practically Applied**



**AneticAid**

# STANDARD FEATURES

**Supremely manoeuvrable, compact, light weight and robust with maximum patient comfort built in: ideal for patient transport, examination and intubation in the OR and any other clinical area.**

With the backrest raised, it is little more than 71" long and it descends lower and is narrower and shorter than nearly all other patient stretchers on the market.



## Trendelenburg tilt

Extremely smooth and controlled mechanism, even at lowest height

## Vertically operated side rails

Store away under mattress for a virtually zero transfer gap. Spring counterbalancing allows effortless ease of operation

## 5th wheel steering

Unique mechanism ensures maximum maneuverability and easily guided by one person

## Patient comfort

Exceptional low height capability and store away side rails allow easy mount/dismount particularly for frail and elderly patients. Wide patient platform and pressure care mattress for patient comfort. Optimum ground clearance for patient hoist access



## Pushing handles

Four ergonomically positioned pushing handles for supreme maneuverability

## Gas assisted backlift

Unique design, precision controlled, easy lift system for enhanced user and patient comfort



## K8 pressure care mattress

Manufactured from temperature-sensitive viscoelastic foam that molds to the patient's body shape, relieving pressure by distributing the load as evenly as possible. Outer material is antimicrobial, water resistant, vapor permeable and latex-free



## Infection control

Smooth molded surfaces and acrylic capped vac formings enclose stretcher mechanics to minimise dirt traps. Anti-bacterial properties and innovative design for easy cleaning

## Fixed transfusion pole

Innovative quick release safety locking mechanism, integrates into either head or foot end of side rails for easy access and storage

## Hydraulic lift mechanism

Exceptional low height capability, adjusted by independently operating pedals on both sides of stretcher



## Casters

Polyurethane tires offer low rolling resistance for high maneuverability, 150mm casters fitted as standard, 200mm optional

## Cylinder trough


Side loading for ease of access, suitable for wide range of cylinder options including F, E or CD

Special bariatric version of this model also available – see page 6






**Minimise risks for practitioners dealing with increasing patient weights with powered positioning: the whole stretcher platform or the backrest can be raised or lowered electronically from the panel at the foot end.**

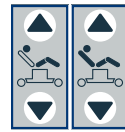

An additional powered feature places the stretcher into the KneeFlex position for easy elevation of the lower limbs for patients having leg and foot procedures – position also provides additional comfort, particularly for patients undergoing abdominal procedures.

There's also a special CPR feature: return the stretcher swiftly to a horizontal, optimal height position for CPR treatment.

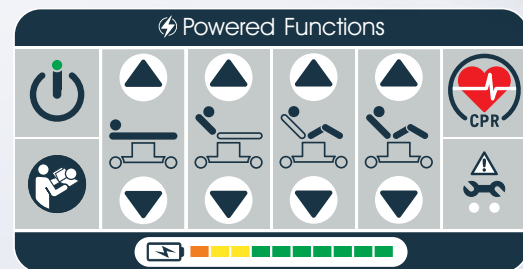

### Electrically powered functions

Easy to operate using the positive feel push button controls on a panel at the foot of the stretcher

Manual override systems as standard


### Platform height – electrically powered raise and lower mechanism

Offers smooth movement to position the stretcher at any height in its range – 18.7" to 30.5" (excluding mattress)




### Backrest – electrically powered raise and lower mechanism

Positions the backrest anywhere between the vertical and horizontal for maximum patient comfort

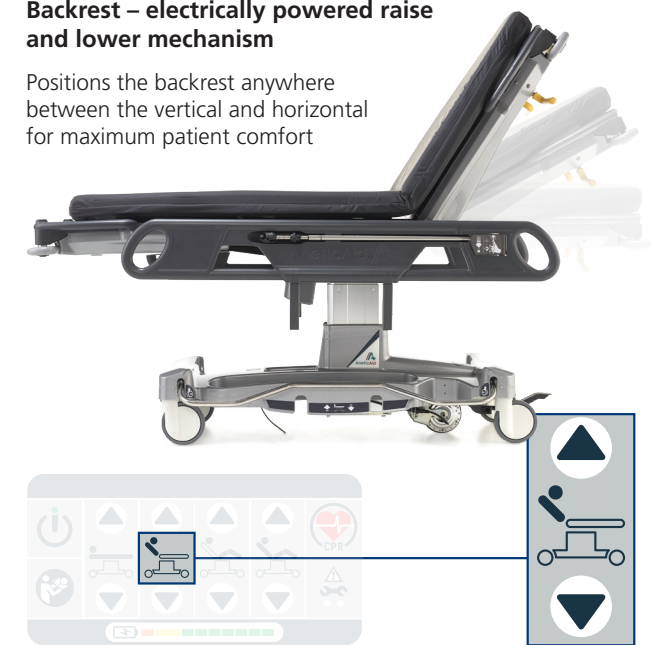

Lights on the power panel illuminate to refer the user to the manual if the equipment is being used incorrectly or to indicate when maintenance is required

**If a patient needs CPR treatment, there is no time to lose: the powered CPR feature returns the stretcher swiftly to a horizontal, optimal height position for CPR treatment**

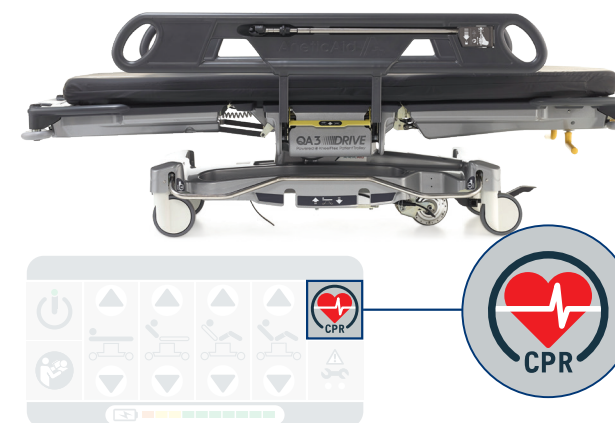

### Trendelenburg – hydraulically powered tilt mechanism

Extremely smooth and controlled movement, even at the stretcher's lowest height



Special features for a stretcher that's more powerful, by design



# STRETCHER MATTRESS

# BUILD OPTIONS

## K8 Pressure Care Mattresses

Pressure care and infection control are top priorities in a stretcher mattress – and Anetic Aid’s independently tested K8 pressure care mattress addresses both of these vital factors.

At its core is temperature-sensitive viscoelastic foam, uniquely specified by Anetic Aid. This means that the mattress can be thinner, while still offering the same, if not better, pressure distribution than a thicker, less sophisticated foam.

This is particularly important with tables and stretchers, as it means the overall height can be kept as low as possible for easy access by frail or elderly patients.

A layer of thermally responsive gel-feel memory foam is then added, further relieving pressure by distributing the patient’s weight as evenly as possible across the surface area – whatever their position.




## Infection control

Both of the foam elements of the mattress have antibacterial and fungistatic properties, while its latex-free, 4-way stretch polyurethane (PU) outer cover is antimicrobial, water resistant and vapor permeable.

This means it is not only contamination-resistant, but it also allows the skin to breathe – for greater patient comfort.

Finally, a unique ultrasonic sealing system effectively welds the seams of the outer cover material together to give an almost unbreakable seal, making the K8 a class leader in infection control.

Description	Code
QA3 K8 Pressure Care Mattress with Zip Cover	21186
QA3 K8 Pressure Care Mattress – Replacement Waterproof Zip Cover	21186-B1

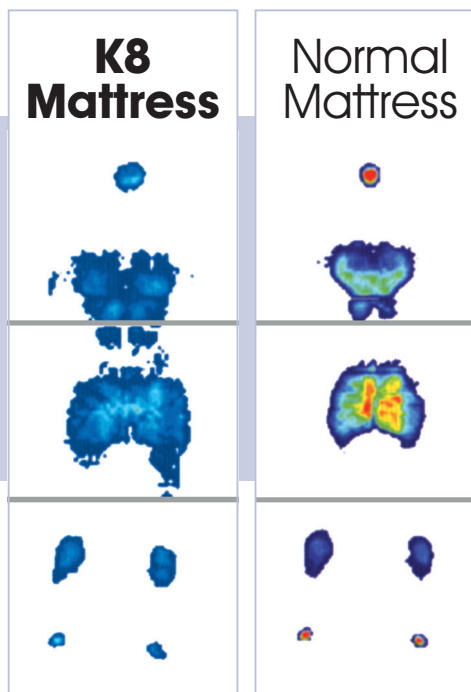

- High risk pressure care mattress
- Independently tested
- Latex-free
- Antibacterial / vapor permeable
- Welded seams
- X-ray translucent



## QA3 <800> for Bariatric Patients

Increased platform / mattress – width: 800mm



Again recognizing the trend for increasing patient weights and sizes, we offer special Bariatric versions of the QA3.

The extra wide (800mm) platform significantly increases the surface area for added safety and comfort for Bariatric patients.

As with all of our stretchers, the QA3 <800> is fitted with the K8 pressure care mattress.



All models of the QA3 have a weight capacity of 705 lbs and can be fitted with this increased platform width

Description	Code
QA3 Build Option <800> 800mm Extended Width Patient Platform	22160

Supports a wide range of surgical specialties



# ACCESSORIES

**QA3™**  
Patient Stretcher

**QA3**  
Powered ⚡ Patient Stretcher



**Monitor Shelf**  
with Removable Refreshment Tray



**Removable Refreshment Tray**



**Push / Pull Bar**



**Foot-end Extension**  
with 4" Pad



**Foot-end Mattress Protector**



**Paper Roll Holder**



**Patient's Storage Container**

SPECIFICATIONS	
QA3 PATIENT STRETCHER (Product code: 21110)	Minimum Height <b>*18.7"</b>
	Maximum Height <b>*30.5"</b>
QA3 POWERED PATIENT STRETCHER (Product code: 21116)	Overall Length <b>85"</b>
	Overall Width (side rails up) <b>33"</b>
	Overall Width (side rails down) <b>29"</b>
	Mattress Depth <b>3"</b>
	Patient Weight Capacity <b>705 lbs</b>

\*NOTE: Measurement from floor to litter top excluding mattress



DESIGN & INNOVATION



MADE IN THE UK



GLOBAL DISTRIBUTION

Innovative **Medical** Technology - **Practically Applied**

**Anetic Aid**, USA + 817 846 5242 [aneticaid.com/us](http://aneticaid.com/us) [steve.smargiasso@aneticaid.com](mailto:steve.smargiasso@aneticaid.com)



**AneticAid**

ISSUE 1 - APR 2021 US