

QA4TM

Mobile Surgery System



The widely adopted single-platform surgery solution for **transport, treatment and recovery.**

AneticAid[®] 

QA4 MOBILE SURGERY SYSTEM

Whether for Day / Ambulatory Surgery or Mainstream Theatre, the QA4 is the perfect combination of transport trolley and operating table.

TRANSPORT ||||| TREATMENT ||||| RECOVERY



Moving and Handling

Zero transfers mean risk is minimised for patient and practitioner alike: from the ward, through anaesthetics, to the surgical procedure, and finally to first and second recovery



Infection Control

A single platform for the entire patient journey, combined with easy-to-clean materials with a range of antibacterial and microbial properties, reduces the chance of cross infection



SOCIAL VALUE

Social Value is at the heart of all we do - engaging our community and giving back



MADE IN THE UK

Designed and produced at our UK manufacturing facility using the latest engineering techniques



CARBON NEUTRAL

We're certified "Carbon Neutral" - meaning we're committed to a better and cleaner environment



10 YEAR WARRANTY

Unique manufacturer-backed Lifetime Warranty Scheme

STANDARD FEATURES

The QA4 Mobile Surgery System is both patient trolley and operating table: a total solution for modern day treatment centres and regular theatres alike.

Removable articulating headrest

Easily swapped for specialist versions such as the ophthalmics headrest



Accessory Side Bar

EU sized accessory side bar measuring at 25 x 10mm



Infection control

Smooth moulded surfaces and acrylic capped vac formings enclose trolley mechanics to minimise dirt traps. Anti-bacterial properties and innovative design for easy cleaning

Vertically operated safety side rails

Store away under the mattress, improving the ease of mounting and dismounting the stretcher for the patient.

Moving and handling

Zero transfers mean risk is minimised for patients and practitioner alike: from the ward, through anaesthetics, to the surgical procedure, and finally to first and second recovery.

Press on / Press off 5th wheel steer

Unique mechanism ensures maximum manoeuvrability and easily guided by one person.

Simultaneous braking low rolling resistance castors

Polyurethane tyres offer low rolling resistance for high manoeuvrability. 150mm diameter central total lock castors fitted as standard, helping to offer a rigid and reliable work platform.



Removable rotating side rails

Providing safety for transport, but removable for maximum surgical access.

Removable adjustable height transfusion / infusion stand

Innovative quick release safety locking mechanism

Removable light weight leg section

Offering maximum surgical access for gynae and urology procedures

O² Cylinder trough

Side loading for ease of access, suitable for wide range of cylinder options including E, F, CD and HX.

Space saving design

With the back rest raised and leg section removed, the system has an incredibly small footprint - making it a comfortable and manoeuvrable recovery or transport position and easier to store when not in use.

QA4™
Mobile Surgery System



Pressure care mattress

Manufactured from temperature-sensitive gel-feel foam that moulds to the patient's body shape, relieving pressure by distributing the load as evenly as possible. Outer material is anti-microbial, water resistant, vapour permeable and latex free. Now 5% resistant to chlorine and other disinfectant compounds.



SPECIFICATIONS		Code: 21300	
Minimum height	*605mm	Mattress width	600mm
Maximum height	*1005mm	Mattress depth	75mm
Overall length	2100mm	Trendelenburg	20°
Overall width	800mm	Reverse Trendelenburg	12°
Castor diameter	150mm	Lateral tilt	±12°
Weight	160kg	Backrest articulation	0 to 65°
Mattress length	2040mm	Headrest articulation	+25 / -30°

*NOTE: Measurement from floor to platform top excluding mattress

A unique combination of manoeuvrability and stability



SOCIAL VALUE



MADE IN THE UK



CARBON NEUTRAL

POWERED FEATURES

Top of the range, the QA4 Powered Mobile Surgery System is ideal for many surgical procedures.

With five remotely-activated powered functions – raise and lower, back rest, traversing top, lateral tilt and Trendelenburg – the trolley offers full 'C' arm coverage and superior surgical access. It also offers exceptional stability with a patient weight capacity of 250kg.

250
kg

250kg Weight Capacity

The weight capacity offers exceptional stability.

Removable remote handset

Easy to operate using the hand control positive feel push button controls.



Patient platform height raise and lower



QA4TM Mobile Surgery System

Lateral tilt + - 12°



Backrest 0 to 65°



Trendelenburg 20° / Reverse Trendelenburg 12°



Head-end / Foot-end Traversing platform



SOCIAL VALUE



MADE IN THE UK



CARBON NEUTRAL

Effortless positioning at the touch of a button

QA4 MOBILE SURGERY SYSTEM

The versatility and superior surgical access of the QA4 Mobile Surgery System make it ideal for a wide range of procedures, including general surgery, gynaecology, urology, endoscopy, ophthalmics, dental and ENT.

ENT (Ear, Nose & Throat)



Head and Neck / OMFS



Ophthalmology



Standard articulating head section or ergonomically designed dual-articulating head section for movement in three different planes

Breast surgery

Can be fitted with Reison arm supports for superior positioning and stability



General surgery

Variable height, Trendelenburg and Reverse Trendelenburg positioning, as well as lateral tilt, make the QA4 ideal for laparoscopic procedures such as cholecystectomy (gallbladder removal)



Urology



Unprecedented surgical access for urologists, plus dedicated Perineal Instrument Tray accessory



General / Orthopaedic / Vascular

Optional leg abduction supports for varicose vein procedures

Gynae/Obstetrics



Supports a wide range of surgical specialisms



SOCIAL VALUE



MADE IN THE UK



CARBON NEUTRAL

QA4 MOBILE SURGERY SYSTEM

Mobile Surgery Systems are increasingly being used for a wide range of orthopaedic procedures: the QA4 can be fitted with a variety of attachments including arm tables, lateral and pelvic supports as well as specialist attachments for knee and shoulder arthroscopy.



Attachments are available for lateral, supine and prone positioning



Knee arthroscopy

Standard leg section removed, making way for specialist leg supports such as the Knee Stabiliser, available in our Operating Table Accessory range



Shoulder surgery

Complete surgical access with shoulder suspension system, or specialist shoulder attachment with back support and headpiece to secure the patient's airway



SURGICAL VERSATILITY

C-arm

The QA4's platform traverses to both head and foot ends, providing comprehensive C-arm access. We also offer a Retrograde Extension accessory for ERCP retrograde tomography



SOCIAL VALUE



MADE IN THE UK



CARBON NEUTRAL

Specialist attachments for orthopaedic surgeries

QA4 MOBILE SURGERY SYSTEM

BUILD OPTIONS

Model specific build option



QA4™ Powered ⚡ Mobile Surgery System

Description	Code
Compact Pushing Handles	22345



Headrest Full Width

Description	Code
Headrest - Full Width Articulating - with Pressure Care Mattress	22322
Headrest - Full Width Articulating - with EU Side Bar and Pressure Care Mattress	22331



Headrest – Dual-articulating

Description	Code
Headrest - Dual-articulating Head Positioner - with Shaped Cushion and Neck Plate Pad	22325

Oxygen Cylinder Bracket

Description	Code
Oxygen Cylinder Bracket - for CD Sized Cylinder	21393



Department Identification

Description	Code
Department Identification - Vinyl Transfer Lettering	21175





Lift-Assist™ Beach Chair

Description	Code
Lift-Assist™ Beach Chair - includes reusable pad set and intubation plate & pad	10786

Customisation

Build Options (which must be specified at the time of ordering) enable you to customise your QA4 with special features – including swapping standard Headrests and Leg Sections for those designed for particular surgical specialisms.



Leg Section – Standard with EU Side Bar

Description	Code
Leg Section - Standard Light Weight - with Full Length EU Side Bar	22358



Leg Section – Articulating with EU Side Bar

Description	Code
Leg Section - Articulating - with 105mm EU Side Bar	22354



Leg Section – Articulating

Description	Code
Leg Section - Articulating	22355

Build options must be specified at the time of ordering



SOCIAL VALUE



MADE IN THE UK



CARBON NEUTRAL

QA4 MOBILE SURGERY SYSTEM

Gardiner Hand Surgery System

Description	Code
The Gardiner Hand Surgery System*	21380

*Requires 2 x Clamps
(Code: 10248-B1 / 10248-B2)



Woodlands Arm & Hand Surgery Table

Description	Code
The Woodlands Arm & Hand Surgery Table	21358



Magna Hand Surgery Table

Description	Code
The Magna Hand Surgery Table*	21384

*Requires 1 x Clamp
(Code: 10248)



Extending Scope

Specialist Headrests and Leg Sections, along with a wide range of other Operating Table Accessories, can be purchased at any time to extend the range of surgical procedures that can be performed on a QA4.



Leg Section – Standard with EU Side Bar

Description	Code
Leg Section - Standard Light Weight - with Full Length EU Side Bar	21357E



Leg Section – Articulating

Description	Code
Leg Section - Articulating	21355



Split Leg Positioner - with EU Side Bar

NOTE: CLAMPS SOLD SEPARATELY

Description	Code
QA4™ Split Leg Positioner - for EU Side Bar	21359



Leg Section – Articulating with Side Bar

Description	Code
Leg Section - Articulating - with 105mm EU Side Bar	21354



Headrest Full Width

Description	Code
Headrest - Full Width Articulating - with Pressure Care Mattress	21322
Headrest - Full Width Articulating - with EU Side Bar and Pressure Care Mattress	21322-EU



Armboard

Description	Code
Armboard - Rotatable (10-387) - with Rotating Clamp - for EU Side Bar	10400-EU
Armboard Cushion / Pad (10-382-H) - Medisoft (Viscoelastic) - Tapered Bottom with Guttered Top Surface - 40mm*	10415-B1

*Required when ordering 10400-EU - Armboard - Rotatable



Headrest – Dual-articulating

Description	Code
Headrest - Dual-articulating Head Positioner - with Shaped Cushion and Neck Plate Pad	21325



SOCIAL VALUE



MADE IN THE UK



CARBON NEUTRAL

Wide range of extra items to compliment the QA4 Mobile Surgery System

QA4 MOBILE SURGERY SYSTEM

Lift-Assist™ Lithotomy Stirrups

Accessory Clamps
supplied separately unless
otherwise stated



Quick reference guide:

STANDARD
Handle type: Twist Grip
Boot type: Standard
Lithotomy range: +84° / -33°
Abduction range: +25° / -9°



HEAVY DUTY
Handle type: Twist Grip
Boot type: Enlarged Standard
Lithotomy range: +84° / -33°
Abduction range: +25° / -9°



FINNED ELITE
Handle type: Squeeze Grip
Boot type: Protective Fin
Lithotomy range: +84° / -33°
Abduction range: +25° / -9°



Description	Code
Lift-Assist™ Lithotomy Stirrups - Standard - includes Boot-pads - Max. Patient Weight @159kg / 25 Stone	10755
Lift-Assist™ Lithotomy Stirrups - Finned Elite - includes Boot-pads - Max. Patient Weight @227kg / 36 Stone	10757
Lift-Assist™ Lithotomy Stirrups - Heavy Duty - includes Boot-pads - Max Patient Weight @327kg / 52 Stone	10758

Retrograde Extension

Description	Code
Urology / Retrograde Pyelography / Gynaecology Extension - with EU Side Bar	21397-EU



Perineal Instrument Tray

Description	Code
Perineal Instrument Tray, Multi-pin	10235



Knee Crutch Stirrups

Description	Code
Simplicity™ Knee Crutch Stirrups - Complete with Boot Pads (Supplied without Clamps) Pair	10765



Oxygen Delivery Bar and Drape Screen

Description	Code
Oxygen Delivery Bar & Drape Screen	21365



Side Rail Covers

Description	Code
Side Rail Cover - PADDED - for Single-piece Side Rail - GREY Material - PAIR	21347
Side Rail Cover - PADDED - for Single-piece Side Rail - CHILD PRINT Material - PAIR	21348



Push / Pull Bar

Description	Code
Push / Pull Bar	21390



Lateral Arm Supports

Description	Code
Lateral / Prone Arm Support c/w Pad	10773



Customise your surgery system with specialist accessories



SOCIAL VALUE




MADE IN THE UK



CARBON NEUTRAL

OPTIONAL EXTRAS

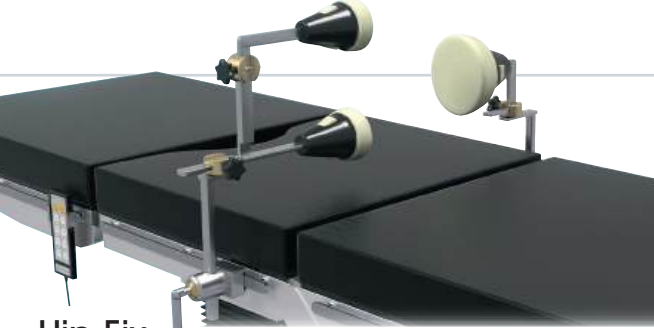

Hip Fix	
Description	Code
Hip Fix™ (2015-2000) 3-point Fixation System - includes sacrum positioner & vertically mounted iliac crest positioner. Single-use pad set & clamps (two required per system) sold separately	10470
Hip Fix™ Rotary Clamp (10-305-P) - Hook-on for UK Side Bar with Flat Pushing Plunge	10470-C1
Hip Fix™ Rotary Clamp (10-305-P) - Hook-on for EU Side Bar with Flat Pushing Plunge	10470-C2
Hip Fix™ Single-use Pad Set - includes one sacrum & two iliac crest positioner pads (Pack of 25)	10470-P1



Clamps	
Description	Code
Clamp - Basic Blade Handwheel - for UK and EU Side Bar	10761
Clamp - Rotary Socket - Stainless Steel - Hook-on - Universal for UK / EU / US Side Bar	10248
Clamp - Universal (Socket & Rail) - Handwheel - for UK / EU / Swiss Side Bar	10250



Lateral Brace 3-Piece

A product by Baxter®, this Lateral Brace comes with extended posts that adjust for use on a variety of patients. It stabilises patients on OR tables, especially for total hip and other procedures requiring optimum lateral stability.

As well as being appropriate for most makes of operating table, the Lateral Brace can be used with the QA4™ Mobile Surgery System - the complete multi-purpose solution for patient transport, treatment and recovery.

Description	Code
Lateral Brace 3-Piece c/w Rectangular x2, Circular x1 Supports	10770

QA4 Mobile Accessory Stand

Accessories need storage when not in use – and we have the perfect solution.



Description	Code
QA4 Mobile Accessory Stand - with EU Side Bar	21370

QA4™ MOBILE SURGERY SYSTEM

The QA4™ isn't just the "hidden standard" of UK day surgery, it's the proven system in enabling safe and efficient surgical work-flow.



MORE CASES



SAFER WORKFLOW



HIGHER EFFICIENCY




Subject to annual servicing



We offer the unique and popular 'Lifetime Warranty Scheme'

Take out a **AnetiCare** service package for regular and predictable maintenance costs that includes a **10 YEAR WARRANTY**. Covers engineer's time and most parts. Contact us to learn more about what an **AnetiCare** service package can offer you.

Available via your Anetic Aid partner – find them at aneticaid.com/distributors.

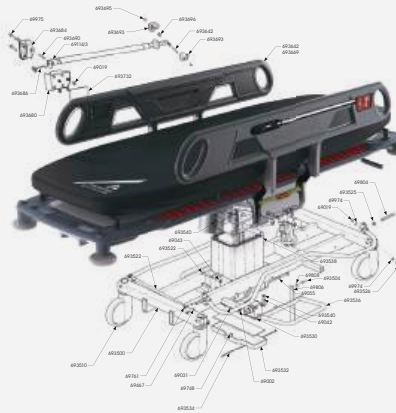

SOCIAL VALUE



MADE IN THE UK



CARBON NEUTRAL

Value for money: minimum downtime, maximum peace of mind



QA3 Ophthalmic Stretcher:

Outstanding surgical access - exceptional patient comfort.

PLUS

our other product areas



Scan this code to see our **QA4's** journey through day surgery!



SEE
ALSO



For more information on any of our products or service contracts, call:
sales@aneticaid.com **+44 (0)1943 878647** aneticaid.com



SOCIAL VALUE



MADE IN THE UK



CARBON NEUTRAL

Innovative **Medical** Technology - **Practically Applied**



Anetic Aid Ltd, Sapper Jordan Rossi Park, Otley Road, Baildon BD17 7AX, UK